


Product range changes as of 1 July 2018

Our product range is always changing. The following pages give a summary of the product changes and novelties of the past half year. The adjusted price lists can once again be accessed by logging in to the website. Have you forgotten your login details? Not a problem. We'll be happy to reissue this information!

Opening of Life & Mobility France

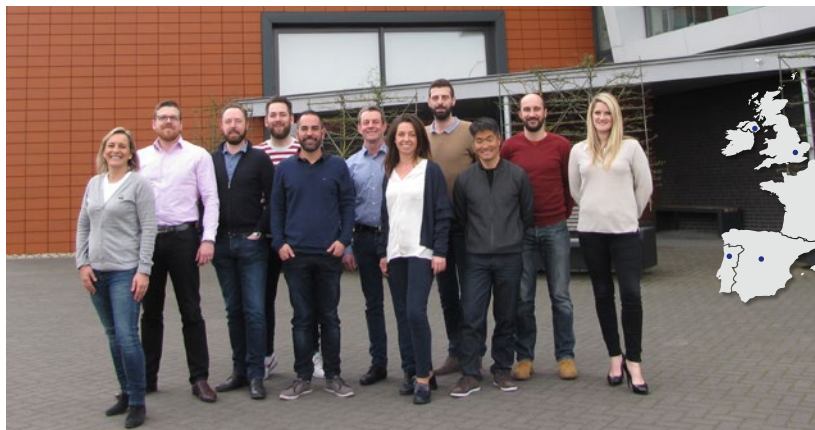


After the opening of our branch in Berlin, our new branch close to Lyon was officially opened last month. This step means that Life & Mobility is now also represented in France. The enthusiastic French team of product advisers recently visited Doetinchem for an update on all the novelties and to exchange experiences with their Dutch colleagues. Stand-up wheelchairs hold an important position in the French product range. A complete overview and more information on the French and German product ranges are available on the Life & Mobility website; www.life-mobility.com.

Online Parts Catalogue

For Life & Mobility, user convenience and growing demand for spare parts were good reasons to further develop the Online Parts Catalogue. For example, you can now easily order a part from the system with your phone. Important improvements are:

- Spare parts for electric wheelchairs have been added
- It is possible to place orders directly from the system
- The system shows prices with the parts
- Suitable for Windows, Android and iOS operating systems
- Suitable for all screen formats, including tablets and smart phones

You can access the Online Parts Catalogue via the link on our website. To order parts, you will need to log in.

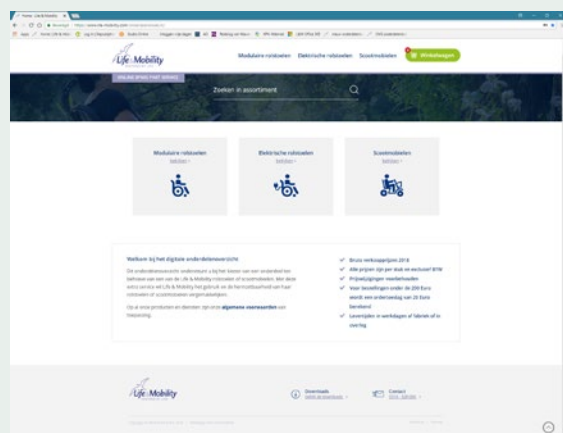

Kerkrade Specialty Centre

Individual modifications are an important part of our service. The Life & Mobility branch in Kerkrade has become a Specialty Centre for powered mobility. Here, our specialists can turn your special requirements for electric wheelchairs into safe solutions. We also have a Specialty Centre for manual mobility. This is at the Life & Mobility head office in Doetinchem. If you have any questions about this, our customer service department will be happy to assist you.



Inspired by life



The perfect Match

The Match is the easily adjustable chassis with which safe combinations can be made with selected seating orthoses. On the basis of user requirements, a number of optimisations have been developed for the Match.

Already implemented

1. The ground clearance below the seating height columns has been increased.
2. The front wheels can be adjusted more easily to be perpendicular to the floor, improving the riding characteristics of the wheelchair.

Upcoming changes

3. The noise made by unlocking the foot brake will be muffled still further, avoiding any potential shock reaction of users or fellow-residents.
4. The bend in the arm rest and its inclusion in the seating plate will be changed in order to create a more durable unit.
5. The information stickers and the explanation of the balance point setting will be clarified so that they contribute towards the correct adjustment of the wheelchair.



Roxx continuous back angle adjustment

Due to changes in the European standards, the mounting of the gas springs of the continuous adjustment of the back-rest angle on the Roxx frame has been improved. The new assembly bracket to be fixed to the frame will replace the single bolt which was previously fixed through the frame. The 'Roxx assembly set for back angle gas spring' can also be ordered separately via the Online Parts Catalogue and can be mounted by the dealer in order to update an existing version.



Our all-rounder: the MyPlus

The MyPlus is our basic wheelchair for adaptations. In the past year, the frame has been redesigned, taking account of the many adaptations already made in the past. The most common adaptations are already listed in the order form, but different requests can also be submitted via the form. If the 'Ocean Blue' colour is not right for you, that is no problem. The MyPlus can be ordered in any RAL colour required, with no additional charge.



Rollit changes

The Rollit has recently undergone a number of changes. As already known, the Rollit now also has the new anti-tip supports and an upgraded stick holder. Like the MyPlus, the Rollit is now supplied with reinforced leg supports and the semi-transparent sides as a standard. Spare parts for the old anti-tip supports, stick holder and leg supports remain available as usual. These can be found in the Online Parts Catalogue.

Gone from the product range

Plus & Rollit spare parts

From 1 July, supplies by RAL colour can be made for as long as stocks last. The standard colour is black.

Matrixx amputation support, old model

Gone from the product range from 1 July.

Spare parts for Solo TS120 and QS, old models

Gone from the product range from 1 July.

Snellius seating element

Gone from the product range from 1 July.

Silver-grey covers for Solo TS120 and QS, old models

Gone from the product range from 1 July.


Cover for the independently adjustable leg support (OIB cover)

Gone from the product range from 1 July.

On track with the @Work auxiliary drive

On the basis of experience in the market, an update has been implemented of the software and hardware (print plate) of the @Work auxiliary drive. This can now be tailored more closely to the customer's wishes. The improvements are:

1. The addition of course correction ('Side slope assist'), which can be turned on and off with the software and can be adjusted according to the weight of the user. With an unexpected sideways move, the wheels make a correction individually, so that the wheelchair continues to ride straight ahead.
 - course correction 'on' is especially for outdoor use on more uneven ground.
 - course correction 'off' is suitable for one-handed propulsion (indoors) for which the user steers with the foot.
2. The motor's *coasting distance* is now more freely adjustable.
 - In other words, the wheels can continue to roll for longer after each push on the hand rims, making it easier for the user to travel long distances.
 - Or the wheels can be set to coast less, for the use of auxiliary drive in small spaces.
3. The sensitivity ('*response level*') of the hand rims is more easily adjustable for each wheelchair user.

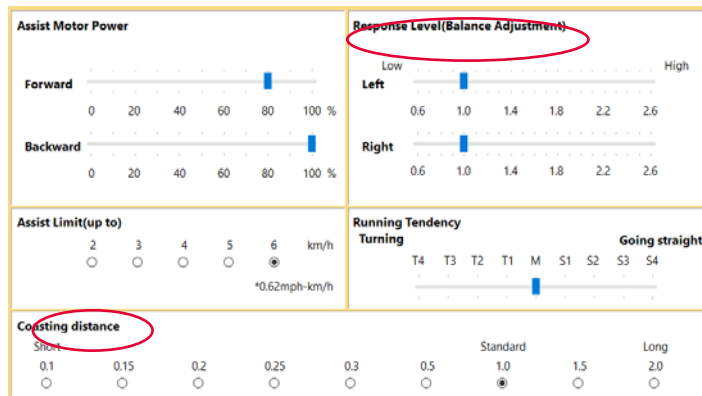

Aura – Manoeuvring lever

For users with little experience with push support, for example, co-residents in an institution or volunteers, it turns out to be difficult in practice to manoeuvre with the motor support of the Aura in smaller rooms or, for example, to position the wheelchair correctly at a table. The Aura push support can therefore optionally be provided with a manoeuvring lever. This option can be used to unlock the brake from the 'off' setting by pressing the lever. This makes it easy to manoeuvre with the wheelchair without having to activate the motors or the wheels getting locked. The Aura equipped with a Manoeuvring lever can be recognised by an extra blue print on the lever and on the control unit.



Mobility scooters VID day

We have implemented a number of optimisations in our mobile scooters range, based on user information:

1. Improved lighting and visibility (Solo)
2. Universal accessories holder (Solo/Mezzo/Primo)
3. Improved steering rubbers (Mezzo/Primo)
4. Improved mirror attachment (Mezzo, with the Primo to follow shortly)

The mobile scooters are now also supplied with a number of extras:

1. Smart on-the-road tips
2. Maintenance booklet
3. The Solo is now supplied with a new robust key with gift box

All modifications will be implemented on the mobile scooters, or can be ordered as accessories, from mid-April!

The above novelties were presented during the VID day for dealers, along with a unique edition of our Solo: the Blue Diamond! This stylish version of the well-known Solo was received with great enthusiasm.





The new Kelvin: take a seat



Our Kelvin seating element, for users who can no longer sit in an ordinary wheelchair, has been upgraded:

1. The seating height of the Kelvin has been lowered, so that users are better able to sit at a table and transfers can be made more easily.
2. In addition, the backrest part of the canvas base is now fitted with a zip which can be used to place new optional support cushions in the desired positions. This allows suitable individual support to be provided while the user is sitting in the chair. The support cushions offer:
 - a. Side support over the entire length of the torso using the long side support cushion
 - b. Local side support using the short side support cushion
 - c. The plastic side plate provides extra stability (is always used in combination with a or b)
 - d. Lumbar support with the lumbar support cushion
3. The backrest angle can be adjusted with an Allen key without unzipping the canvas. This makes it easier to set adjustments.
4. The new, optional back corner cushion can be placed in the lowest corner of the canvas base. This provides extra support for users who require support with a larger backrest angle and is available in 3 widths.
5. The new Kelvin is immediately recognisable thanks to the round zip in the seating element's backrest, combined with the silver print of the Life & Mobility bird.

The Roxx connected!

Together with business software developer SAP, we have started a research project in which the seat and back-rest of the Roxx have been equipped with a number of sensors that we have connected to the Internet. Through this innovation, we can track the sitting behaviour of the wheelchair user. We use this information in the further development of our Matrixx seating system. Together with the Locomotion research agency and the Sensire Den Ooiman care institution, we have also conducted a scientific study of the relationship between sitting (properly) and welfare among wheelchair users. An interesting film has been produced about this and is available on our website. You can view this via this link:
www.life-mobility.com/nl/nieuws/smart-wheelchair-prototype/



Morgan tyres

From now on, the Morgan will be supplied with black tyres as a standard. The black non-marking air tyres fit well with the black frame and so give the Morgan an attractive modern image. Naturally, the grey tyres can still be ordered as an option and massive PU tyres are also available for both colours. The same will soon be available for the Leon. We will keep you informed!



Privacy law

We know a great deal about our customers, but do everything possible to protect that information. If you wish to know how we treat your data at Life & Mobility and how privacy is regulated on our website,

please read our updated statement on the website:
www.life-mobility.com/nl/privacyverklaring.